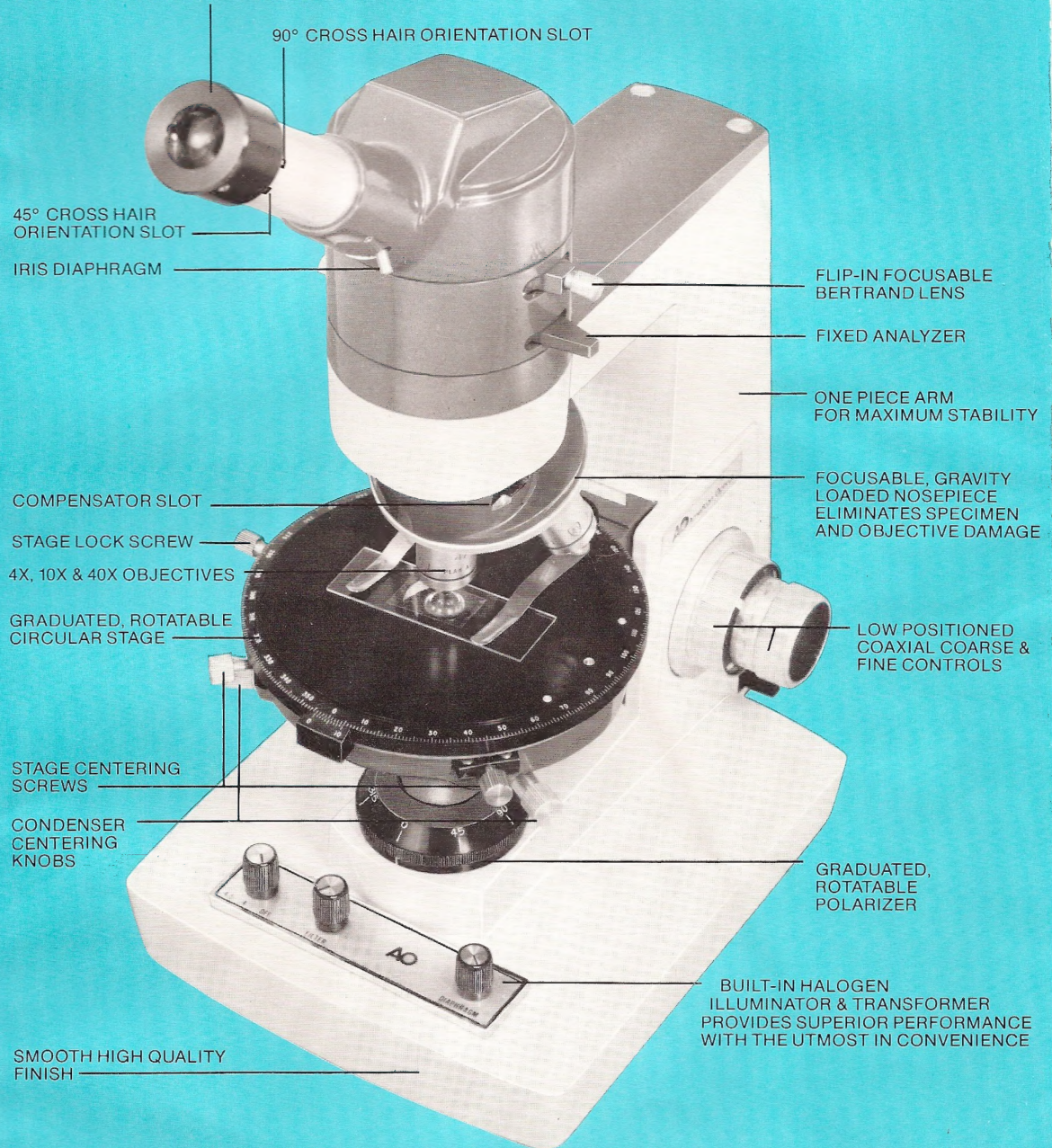




AO POL STAR

POLARIZING MICROSCOPE SERIES 110

10X WIDE FIELD EYEPIECE WITH CROSS HAIR RETICLE
AND EXTRA LARGE 20mm FIELD OF VIEW

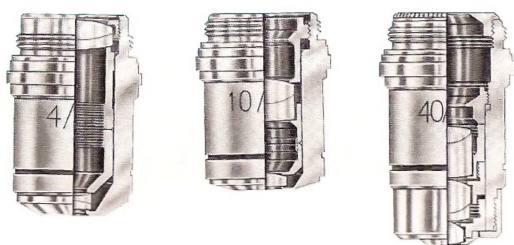


COMPENSATORS (Standard DIN Size – 6.0mm x 20mm)	1886	1/4 Wave Mica Plate
	1887	Full Wave
CONDENSERS & CONDENSER MOUNTS	1091	Auxiliary Condenser
	1092	Condenser Mount
	1087S	Abbe Condenser N.A. 1.25 Strain Free
ILLUMINATORS & LAMPS	1130E	In Base Transformer & Tungsten Halogen Illuminator for 220V, 50Hz Lamp
	1120	Replacement Lamp for 1130 Illuminator 6V, 20W Tungsten Halogen

AO[®] American Optical
INSTRUMENT DIVISION
BOX 123, BUFFALO, NY 14240

- **Extra Wide Field**
- **Cross Hair Orientation**

The AO POL STAR Polarizing Microscope is equipped with a wide field 10X focusable eyepiece offering a full 20mm field of view. It also has cross hairs which can be oriented in either a 45° or 90° viewing position for mineral identification.



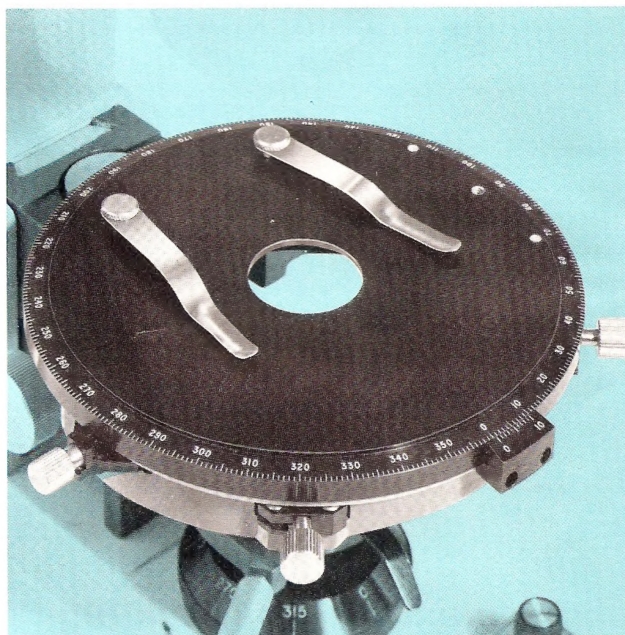
- **Superior Image Quality**
- **Bright, Contrasty Interference Figures**

The performance of any microscope is only as good as its optical system. Therefore, the AO POL STAR Polarizing Microscope comes equipped with a 4X Plan Achromatic, 10X strain free Plan Achromatic, and a 40X strain free Achromatic with N.A. 0.85.

All objectives are parcentered and parfocalized to minimize centering and focusing time. To insure stability and smooth operation, the objective nosepiece is of ball bearing construction. Traditional color coding assures quick positioning of the desired objective into the optical path.

All essential optical surfaces are AMERICOTE™, an AO exclusive coating process which decreases internal reflections, improves light transmission and image contrast.

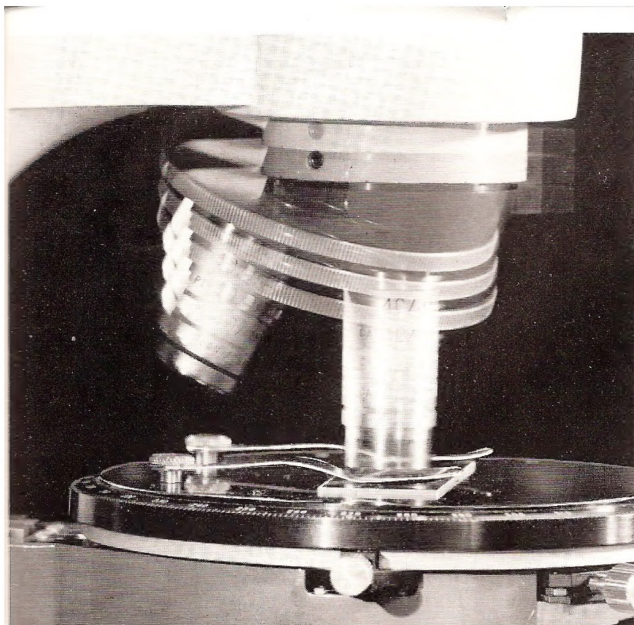
Of utmost importance is the contribution of infinity correction to optical performance. Constant tube length between objective and body is no longer necessary. This exclusive AO design results in a new level of exceptional performance and optical excellence.



- **Graduated**
- **Rotatable**
- **Circular Stage**

The AO POL STAR Polarizing Microscope's 360° rotatable, circular stage allows the specimen to be completely rotated and still remain centered to the optical axis.

The stage is marked in degrees with a vernier scale reading to one-tenth of a degree and is precentered at the factory. Two centering knobs are provided for future centering. A locking screw is also included so the stage can be set in a fixed, non-rotatable position.



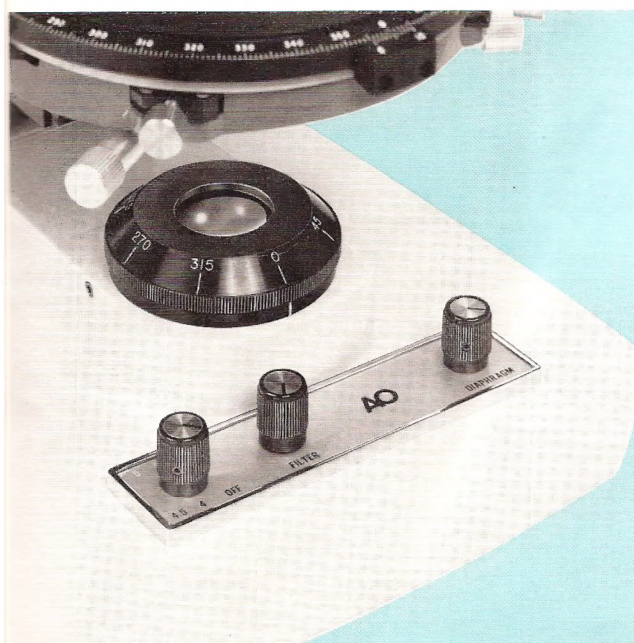
• Unique Focusing Nosepiece Concept

Only the nosepiece assembly moves to focus objective to specimen. Since the stage is not part of the focusing procedure, it is rigidly attached to the stand. The result is a stable, vibrationless system, free from the effects of hand pressure on the stage.

The maintenance-free focusing mechanism is enclosed in the microscope stand assembly, protected from dirt or damage. Coaxial coarse and fine adjustment controls are mounted low on the stand for comfortable, convenient operation.

The nosepiece is gravity loaded. There is no positive downward force. If the objective is accidentally lowered onto the slide, the nosepiece assembly can "give," avoiding objective or specimen damage. Movement of the nosepiece on a ball bearing slideway insures smooth positive action, free from backlash.

ES

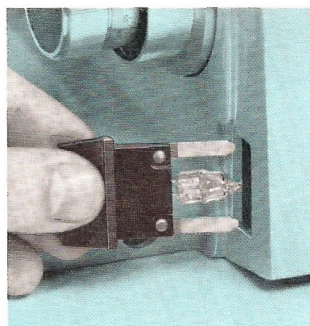


• More Light-Increased Lamp Life • Easy Change Bulb • Clean Burning, High Color Temperature Tungsten Halogen Illuminator

AO's built-in Illuminator is one of the most advanced light sources ever offered on a laboratory microscope. It combines superior performance with convenient operation. All of the following features are incorporated:

- Quick change High Intensity 20W Tungsten Halogen Lamp
- Continuously variable intensity built-in transformer
- Solid State construction
- Built-in field diaphragm and 10% neutral density filter
- Highly corrected lens system

Intensity, Field Diaphragm and Neutral Density Filter controls are conveniently located on the front of the base. Quick, efficient set up of desirable Koehler illumination is assured.



AO
American Optical

The AO Series 110 POL STAR Polarizing Microscope offers transmitted light examinations, including orthoscopic and conoscopic investigations, for industrial, teaching and laboratory applications in:

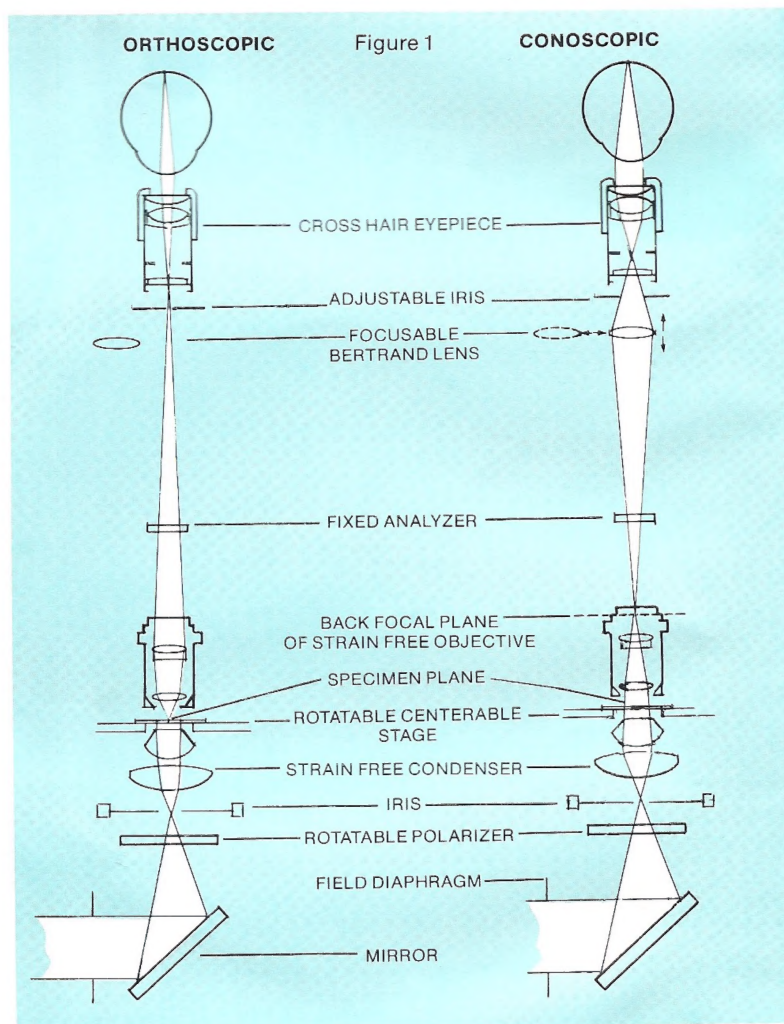
CHEMISTRY • BIOLOGY • GEOLOGY • PETROLOGY • MINERALOGY
METALLURGY • TOXICOLOGY • CRYSTALLOGRAPHY • PULP/PAPER • PAINT
PLASTICS • PHARMACEUTICALS • GLASS • CERAMICS • FORENSIC SCIENCES
POLLUTION CONTROL • GOUT/PSEUDO GOUT



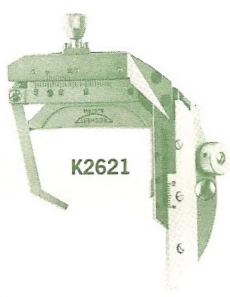
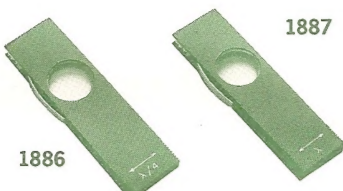
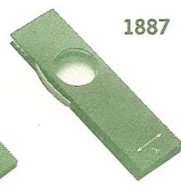



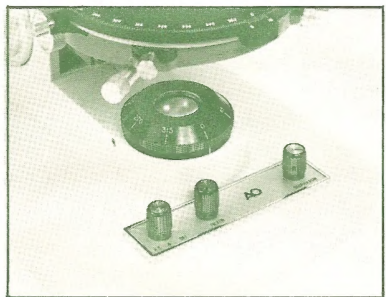
The broad range of information that can be quickly obtained by observing the optical characteristics of a specimen, such as, refractive index, the degree of double refraction, intrinsic color, shape and cleavage, with the AO POL STAR Polarizing Microscope can mean significant time, money and materials savings ... particularly in research and process control areas. It is also ideal as a teaching microscope for courses in Chemistry, Geology, Petrology, Mineralogy, etc.

Figure 1 schematically illustrates the two general systems used to observe, classify and identify substances.

The orthoscopic system is applied to view the sample under crossed polarizer and analyzer. If the specimen is invisible or shows no birefringent properties, it is classified as isotropic; if it exhibits birefringence or is alternately bright and dark or colored upon stage rotation, it is classified as anisotropic.

The conoscopic system is the means for viewing the back aperture of a high power objective and reveals the interference figure produced by an anisotropic specimen. This method permits the classification of anisotropic materials as uniaxial or biaxial; and by insertion of compensators into the compensator slot, determines the nature of birefringence and optic sign.



		CATALOG NUMBER	DESCRIPTION	PRICE
		1882		
		1619		
		K2621		
		1886 1887		
		1881		
		1091S		
		1092		
		1087S		
		1130		
BODY				
1889	Monocular with focusable Bertrand lens and fixed analyzer for POL STAR, only			\$525.00
STAGE				
1882	Graduated, circular, complete with bracket			750.00
K2621	Mechanical Stage Attachment for 28mm x 48mm slides			300.00
CABINET				
1619	Deluxe Cabinet			75.00
COMPENSATORS and POLARIZER				
1886	1/4 Wave Mica Plate			72.00
1887	Full Wave Plate			72.00
1881	Rotatable Polarizer			160.00
CONDENSERS AND CONDENSER MOUNT				
1091S	Auxiliary Condenser			35.00
1092	Condenser Mount			35.00
1087S	Abbe Condenser N.A. 1.25 Strain-free			55.00
ILLUMINATORS				
1130	In-Base Transformer and Tungsten Halogen Illuminator			290.00
1130E	Same as 1130, but for 220V, 50 Hz.			325.00
LAMP				
1120	Replacement Lamp 6V, 20W Tungsten Halogen			7.00

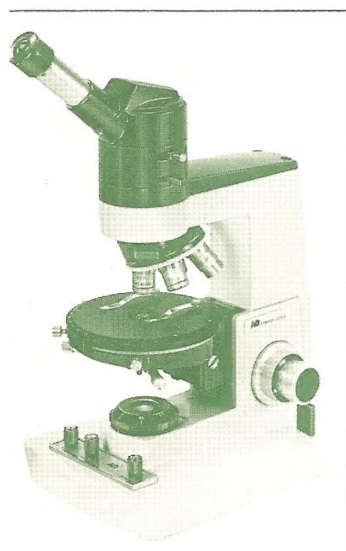
All prices and terms of sale are subject to change without notice.
Applicable taxes are in addition to prices stated.
The equipment supplied may not agree in all details with our description or
illustrations because instruments are subject to modification and improvement.

AO[®] American Optical
INSTRUMENT DIVISION
BOX 123, BUFFALO, NY 14240



PRICE LIST


Effective February 1, 1981

Series H110 POL STAR® Microscope









H110M-POL

STAND & ILLUMINATOR	BODY
 <p>H110</p> <p>1880 Stand, with Quadropole Nosepiece</p> <p>1130 In-Base Transformer and Illuminator, 115V, 60 Hz (220V, 50 Hz available)</p>	 <p>M</p> <p>1889 Monocular</p>
CATALOG NUMBER	PRICE
H110M-POL	\$3100.00

STAGE, OBJECTIVES, CONDENSERS	
 <p>POL</p>	
1882	Circular, Graduated, Rotatable Stage
1890	10X W.F. Eyepiece
1017S	4X Objective
1019S	10X Objective
1891	40X Objective
1091S	Auxiliary Condenser
1092	Condenser Mount
1087S	Abbe Condenser N.A. 1.25, Strain-free

ALTERNATE PARTS AND ACCESSORIES

CATALOG NUMBER		DESCRIPTION	PRICE
 <p>1017S</p>  <p>1019S</p>  <p>1022S</p>  <p>1891</p>  <p>1029S</p>			
 <p>1890</p>			
OBJECTIVES			
1017S	4X, N.A. .12 Strain-free, Planachromatic		\$ 95.00
1019S	10X, N.A. .25 Strain-free, Planachromatic		95.00
1022S	20X, N.A. .50 Strain-free, Planachromatic		150.00
1891	40X, N.A. .85 Strain-free, Achromatic		275.00
1029S	100X, N.A. 1.25 Strain-free, Planachromatic		340.00
EYEPIECE			
1890	10X W.F. with cross-lines		95.00

AO POL STARTM

POLARIZING MICROSCOPE
Series 110



AO[®]
American Optical